


SHOR-LINE®

Inox Series Type 3 Stainless Steel Kennels

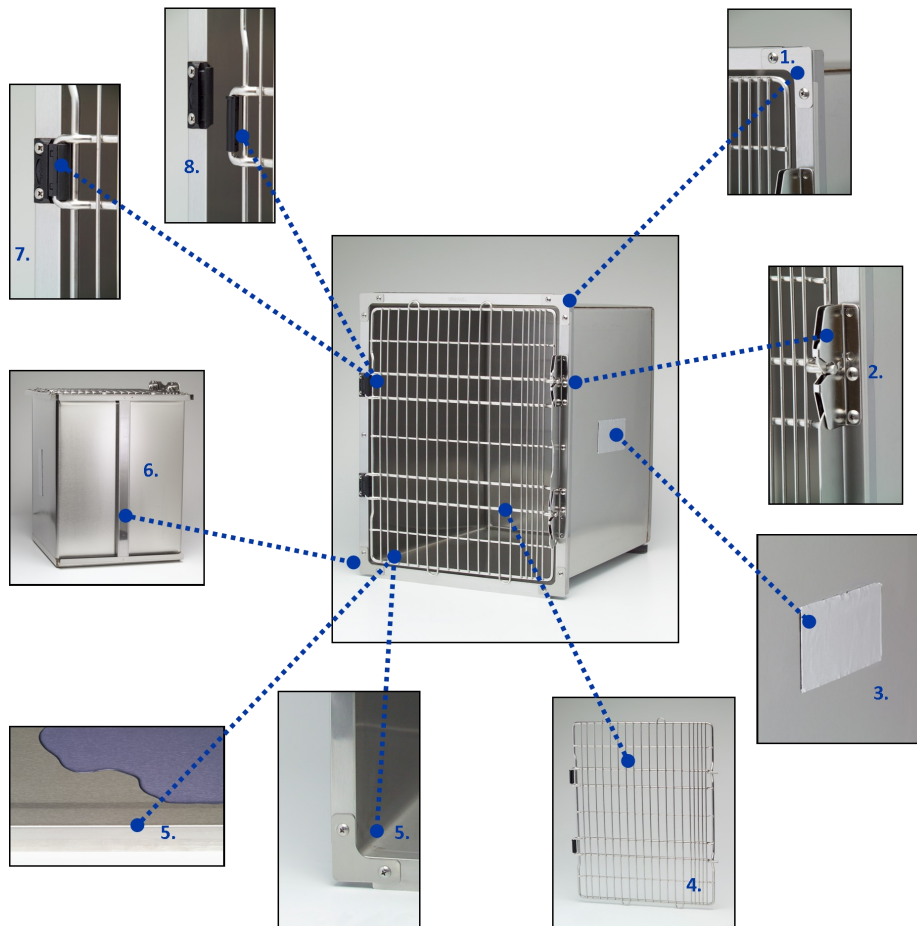

1. Hardware

All Trim Plates and Trim Strips are made from stainless steel with rounded corners. No sharp edges or points.

2. Patented Latch

The latch is designed to be durable, secure and quiet. Acetal bushings provide smooth operation and sound dampening qualities. The concave latch bar design enables the kennel door to be opened with the forearm when handling sedated or awkward animals.

3. Sound Deadening Pads

New sound and vibration absorbing technology which comes standard on all Shor-Line kennels.

4. Kennel Door

The one piece frame and wrap around hinge/latch bars are formed from 9.5mm diameter solid stainless steel rods. All additional horizontal rods are 6mm in diameter and all vertical rods are 4.8mm in diameter. Every intersection of horizontal and vertical rods are welded, providing incredible strength.

5. Fully Rounded Corners and Fluid Retaining Lip

No 90° corners to trap germs and debris. A smooth rounded cove transitions the back and all internal sides of the kennel. A slightly raised, fluid retaining lip prevents drainage out of the front of the kennel.

6. Stainless Steel Frame

All frame members are made from heavy duty stainless steel to provide you many years of solid structural support from your kennel assemblies.

7-8. New Hinge Design

Our new hinge has been designed and developed for strength and utility. The hinge is constructed of 14 gauge stainless steel that is encapsulated by an injection moulded, fibreglass reinforced thermoplastic polymer.

The snap together bushing on the door frame is also moulded from thermoplastic polymer with permanent, self lubricating properties resulting in a silky smooth and silent swing.

The new hinge design also allows for the quick and easy removal of the door for cleaning and reattachment without tools.

Our doors, latches and hinges are completely transposable which allows you to change your kennel from a right latch to a left hand latch, or vice versa at any time.

Quality Equipment

An Investment for Life!